

FOR TRUSTWORTHY BP MONITORING AT HOME

TRUSTCHECKTM BPM^{2.0}

DIGITAL UPPER ARM BLOOD PRESSURE MONITOR

3 YEARS*
EXTENDED
WARRANTY



- WHO INDICATOR
- UNIVERSAL CUFF
- BLUE LED LARGE DISPLAY
- 99 TESTS MEMORY EACH FOR 2 USERS
- IRREGULAR HEARTBEAT DETECTION
- FACILITY FOR TYPE-C USB POWER CABLE

UNIVERSAL CUFF (Fits 22-42cm arm)




REGULAR CUFF (Fits 22-32 cm arm)



TRUSTCHECKTM BPM^{2.0}

DIGITAL UPPER ARM BLOOD PRESSURE MONITOR

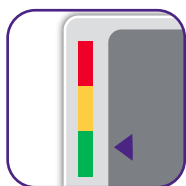

WHO INDICATOR

Makes it easy to understand results



UNIVERSAL CUFF

Fits multiple arm sizes 32-42CM




BLUE LED LARGE DISPLAY

for easy of viewing results



99 TEST MEMORY

each for 2 USERS to store readings of every user separately

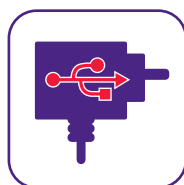

FACILITY FOR TYPE-C POWER CABLE

To use the meter without batteries




IRREGULAR HEARTBEAT DETECTION

Can help in detection of pulse irregularities at an early stage.


TECHNICAL SPECIFICATION

Model	TRUSTCHECK BPM 2.0
Type	Upper Arm Digital Blood Pressure Meter
Display	Digital LCD Display, 58 x 67 mm
Measurement Range	Pressure 0-299 mmHg Pulse 40-180 beats / min
Accuracy	Static pressure: ± 3 mmHg Pulse: ± 5 % of the reading
Memory	99 Results x 2 Users
Dimensions	Approx. 129(L) X 100(W) X 68(H) mm
Weight	230 gm approx. (Excluding batteries)
Cuff Circumference	22-42 cm
Cuff Pressure	0-299 mmHg
Power Source	Internal-DC 6V (4*AA) / External-DC USB 5V
Operating Conditions	Temperature: 10°C to 40°C Humidity: 15% to 90%RH
Storage And Shipping Conditions	Temperature: -20°C to 55°C Humidity: $\leq 95\%$ RH
Measurement Method	Oscillometric Method

WHO INDICATOR INTERPRETATION

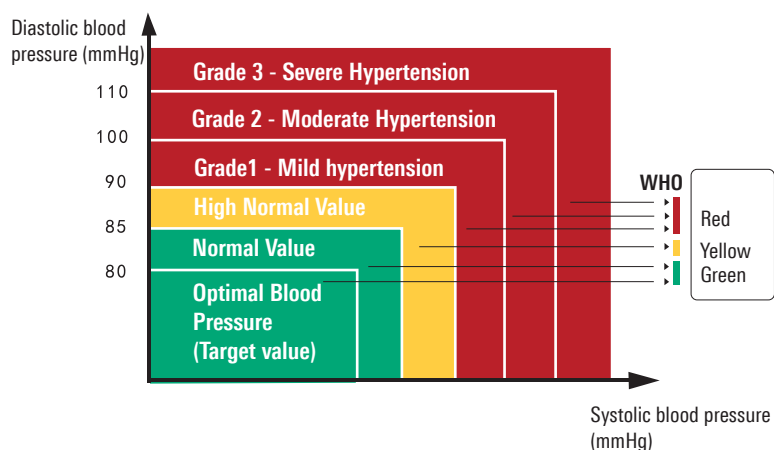

For more information contact:

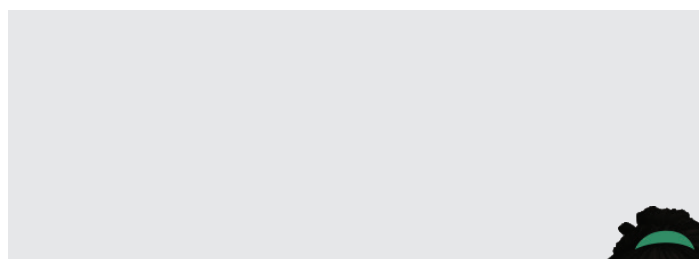

www.arkray.co.in



7777053311



ARKRAY Healthcare Pvt. Ltd. 701-702, Opulence, 6th Road, TPS III, Santacruz (E), Mumbai - 400055, Maharashtra, India.